




VersaSpeaker™

Model VSP-1
Manual: English

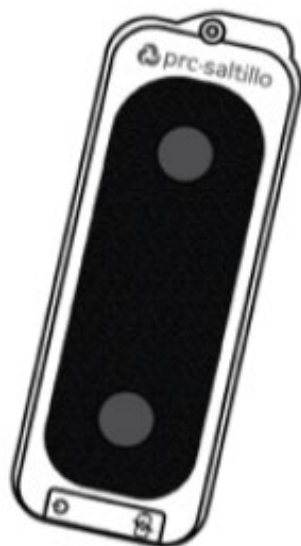

DOC100321v101

Copyright and Trademark Information

© 2025 PRC-Salttillo. Non-commercial reprint rights for clinical or personal use granted with inclusion of copyright notice. Commercial use prohibited; may not be used for resale. Contact PRC-Salttillo for questions regarding permissible uses.

AAC Coaching Corner™, AAC Language Lab®, AAC Learning Journey™, Accent®, AeroSwitch®, AeroText™, ChatEditor™, ChatFusion™, ChatPoint™, ChatSMS®, ChatWrap®, CoreScanner™, Dialogue®, Empower®, Engage®, Essence®, Freedom Switch™, iSharePRC®, LAMP™, LAMP Words for Life®, Look®, Minspeak®, NovaChat®, NuEye®, NuPoint®, NuVoice®, PASS™, PRiO®, Realize Language™, Semantic Compaction®, Timocco-on-Accent™, TouchChat®, UNIDAD®, Unity®, VersaChat™, VersaEye™, VersaSpeaker™, VersaWrap™, Via®, VocabPC™, and Xtreme™ are trademarks of PRC-Salttillo.

AirDrop®, App Store®, Apple®, Apple Books®, FaceTime®, Guided Access®, iCloud®, iPad®, iPad Air®, iPad mini®, iPadOS®, iPad Pro®, iPhone®, iPod®, iTunes®, Lightning®, Mac®, Macintosh®, Safari®, Siri®, TouchID®, and TrueDepth® are trademarks of Apple Inc., registered in the U.S. and other countries and regions.

The *Bluetooth®* word mark and logos are registered trademarks owned by the Bluetooth SIG, Inc. and any use of such marks by PRC-Salttillo is under license.

iOS® is a trademark or registered trademark of Cisco in the U.S. and other countries and is used under license.

Wi-Fi® is a registered trademark of Wi-Fi Alliance.

All other trademarks and copyrights are the property of their respective owners.

Versa Support

Email: versasupport@prc-salttillo.com

Web: <https://www.prc-salttillo.com/versa>

Document Version

DOC100321v100 · 05/04/2026

Table of Contents

Safety Guidelines.....	3
Turn Power On/Off.....	5
Connect to the Device.....	5
Charge the VersaSpeaker.....	6
Light Indicators Reference.....	6
Carrying/Mounting Options.....	7
Option 1: Attach the Lanyard to the VersaSpeaker.....	7
Option 2: Attach the VersaSpeaker to a VersaWrap.....	7
Specifications.....	9
Compliance Information.....	10
International Contact Information.....	12

Safety Guidelines

Cleaning, Maintaining, and Servicing

- ⚠ Do not attempt to disassemble, service, or repair the speaker. If the speaker, charging cable, or lanyard are damaged, call Service for assistance.
- ⚠ Use a microfiber (lint-free) cleaning cloth to wipe the speaker. Never use abrasives or cleaning cloths that could damage the grille or scratch the case.
- ⚠ If the case needs more than a dry wipe-down, dampen a lint-free cloth. “Damp” means wrung out, never dripping. Dry the unit thoroughly with a dry lint-free cloth.

Attaching the VersaSpeaker to a VersaWrap™

- ⚠ Be sure to follow the instructions in this manual if you attach the VersaSpeaker to a VersaWrap. Failure to install the mounting system according to the instructions may result in an injury to the user.

Cables, Straps, or Lanyards

- ⚠ Cables, straps, and lanyards are potential entanglement or strangulation hazards. Please consider this prior to placing these items with users. Never leave small children unattended with the charging cable.
- ⚠ On a regular basis, visually inspect the charging cable and USB-C port for damage or wear. For example, a frayed cable could cause the battery charger/power adapter to short-circuit, posing a hazard to the user.
- ⚠ When unplugging the cable, grasp the connector instead of the cable itself to avoid damage or wear.
- ⚠ If the charging cable becomes damaged, contact PRC-Salttillo Service for a replacement.

Charging

- ⚠ Disconnect the charger when it is not in use to prevent overheating or a fire hazard.
- ⚠ While the product is charging, keep it away from flammable liquids or combustible materials.
- ⚠ Disconnect the charger when it is not in use to prevent overheating or a fire hazard.

Storage

- ⚠ Store the product away from anything that could ignite, including flammable liquids, combustible materials, and direct sunlight.

Interference

- ⚠ Keep the product at least six inches from magnetic-sensitive implantable medical devices such as a pacemaker, implantable cardioverter defibrillator, vagus nerve stimulator, shunt, or stent.
- ⚠ When operating this product in a medical environment, do not use it with any product that is not medically approved. Follow all rules for appropriate wireless product use.
- ⚠ If you have any questions about using your wireless mobile device with an implantable medical device, consult your health care provider.

Water

- ⚠ This product is intended for use in normal situations. It is water-resistant, not waterproof. Do not immerse it in water. Do not spray or splash liquid directly on it and use extreme caution when using it around water. Never use the product when you are actually in the water (for example, a pool or bath).

Heat

- ⚠ Do not place the product in an appliance such as a microwave or oven; near an open flame; on or near a source of heat such as a stove, heater, or furnace vent.
- ⚠ Do not leave the product in a location exposed to direct sunlight or in a high-temperature environment such as a hot vehicle for a prolonged period.

Damage

- ⚠ Do not use the product if it is cracked or broken.
- ⚠ Do not leave the product on an unstable surface from which it could fall and injure the user.

Environment

- ⚠ A hot environment can limit or stop batteries' ability to charge.
- ⚠ Dispose of all electronic components according to local regulations.


Volume/Hearing

- ⚠ Keep volume levels low to prevent hearing damage. Prolonged exposure to high volumes can cause permanent hearing damage. A person can become desensitized over time to high volume levels which could be impairing hearing
- ⚠ Prolonged exposure to loud sounds can lead to tinnitus. Stop using the product if the user experiences ringing or buzzing in the ears.
- ⚠ To set volume to a safe level, adjust the volume to a low setting and then slowly increase the sound until the user can hear it comfortably and clearly.

Turn Power On/Off

To turn power on, press the I/O button.

I/O button

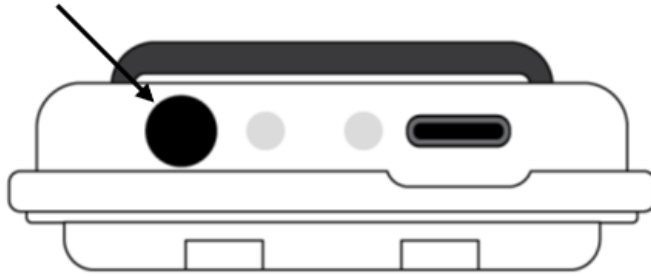

- You will hear a tone that indicates the speaker is powered on.
- The speaker will automatically go into pairing mode.


To turn power off, press and hold the I/O button for 7 seconds.

To perform a hard reset when necessary, press and hold the I/O button for more than 10 seconds.


Connect to the Device

To pair the VersaSpeaker with an iPad, hold down the I/O button for 3 seconds.

The flashing blue light indicates the speaker is in pairing mode.

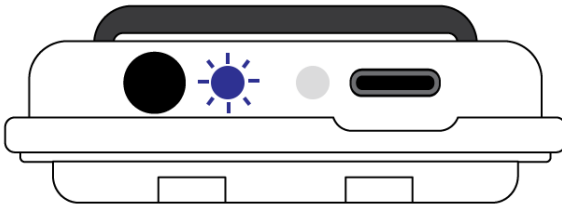


The solid blue light indicates the speaker has successfully paired.

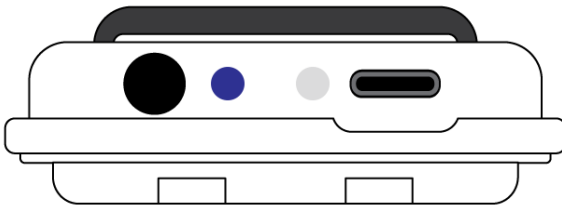

Charge the VersaSpeaker

Plug the charging cable into the USB-C charger port.


- The indicator light should be flashing orange to indicate it is charging.
- The solid orange light indicates the battery is charged (when the charger is attached).

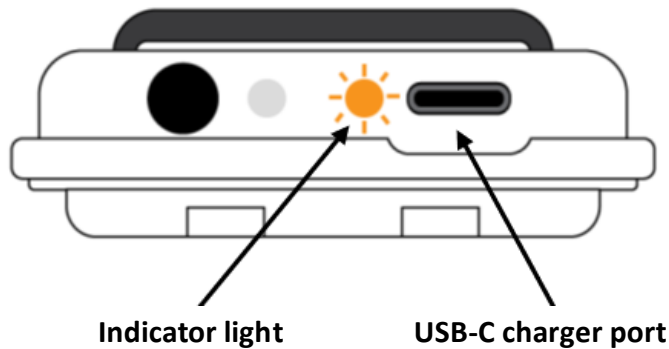

Low Battery: When the orange light flashes twice every 30 seconds, the battery is below 15%. Typical battery life is 9 hours.

- ⚠ While the product is charging, touching the surface of the product with broken skin may aggravate a wound.
- ⚠ While the product is charging, infants or high-risk groups should not touch the surface of the product if there is a chance of burning the skin.
- ⚠ While the product is charging, keep it away from flammable liquids or combustible materials.
- ⚠ Disconnect the charger when it is not in use to prevent overheating or a fire hazard.

Light Indicators Reference

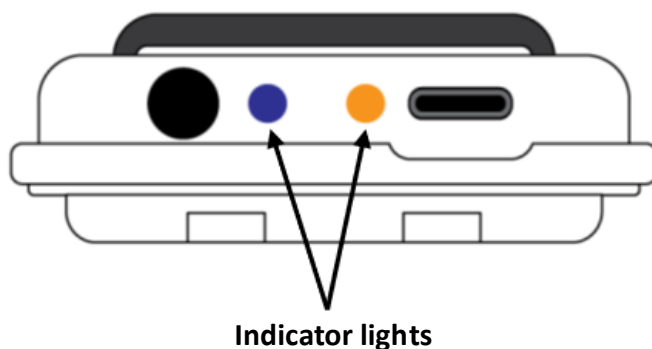

Light	Status
Flashing blue	Pair mode
Solid blue	Paired
Flashing orange	Charging (while plugged in)
Solid orange	Charged
Intermittent flashing orange	Battery level below 15%

Carrying/Mounting Options

Option 1: Attach the Lanyard to the VersaSpeaker

Attach the lanyard that was provided with the VersaSpeaker and wear it to transport the speaker wherever you go.




Option 2: Attach the VersaSpeaker to a VersaWrap

If you purchased a VersaWrap™, you can attach the VersaSpeaker to the wrap and mount the speaker on an iPad®.


Step 1: Remove the Cover Plate

1. Grasp the tab on the right side of the cover plate.
2. Lift the plate to remove it.

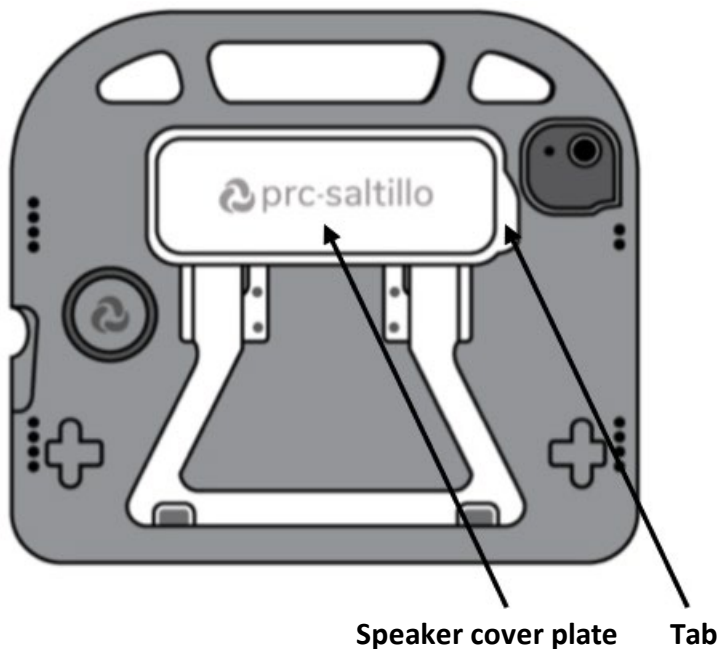


Step 2: Attach the Speaker

1. With the cover removed, align the tabs on the left side and slide the speaker into place.
2. On the right side, press down until the speaker clicks into place.

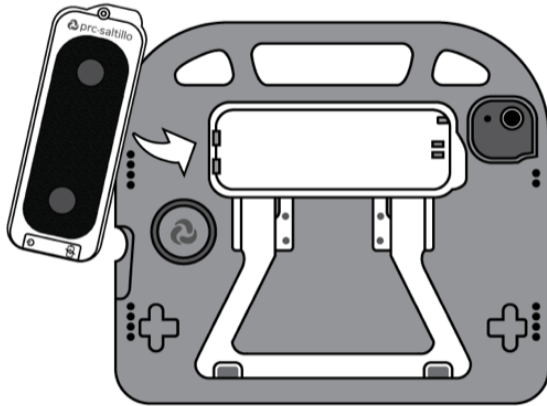

3. On the right side, secure the speaker with the screw that was provided with the speaker.

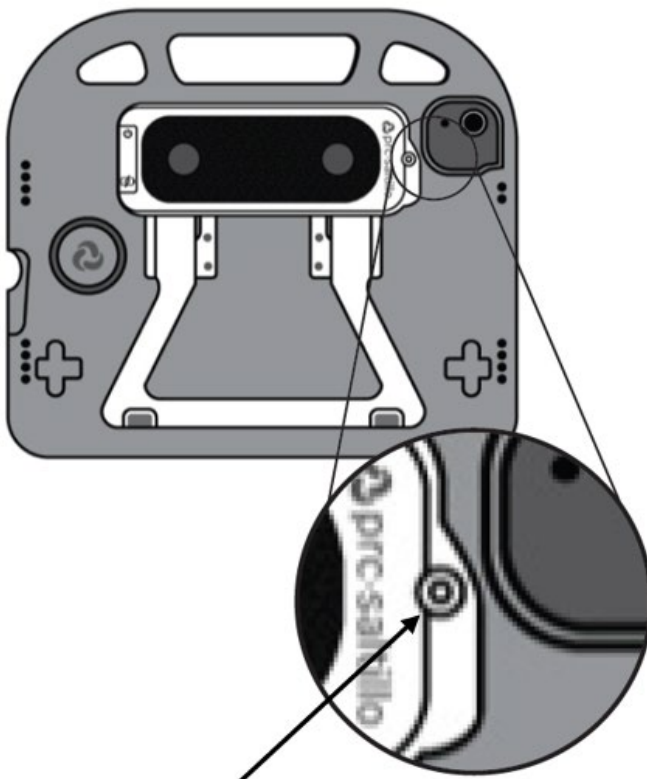

Screw location

Specifications

Weight

9.1 oz.

Dimensions

0.75 in. D x 6.2 in. W x 2.2 in. H

Wireless communication technology/connectivity protocol

Bluetooth

Frequency range

2402 MHz to 2480 MHz

Audio driver type

Dynamic driver

Batteries required

1 chargeable battery included / DC 3.7V

Battery life

8 hours

Ports

USB-C input 5.0V — 2.0A

Power supply/adaptor

Input: 100-240V~50/60Hz 0.4A Max

Output: 5.0V — 2.0A, 10W

Protection against harmful ingress of water

IP44

Ambient temperature

5°C to 40°C (41°F to 104°F)

Included in the package

VersaSpeaker, charging cable, lanyard

Compliance Information

Federal Communications Commission (FCC) Warning

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications to this device not expressly approved by PRC-Salttillo could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radiocommunications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Innovation, Science and Economic Development Canada (ISED) Warning

This device complies with ISED Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

Le présent appareil est conforme aux CNR d'ISDE Canada applicables aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes : (1) l'appareil ne doit pas produire de brouillage, et (2) l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement.

Product Intended Purpose

The product is intended to be used with smartphones, tablets, laptops, or MP3 players indoors, outdoors, and on-the-go.

Contraindications: Use the product with caution if the user suffers from any hearing disorder. Keep volume set to the lowest value at which the user can clearly hear the broadcast sound. Stop using the product if it impairs hearing. For details, see "Volume/Hearing" under "Safety Guidelines" in this manual.

Security

Turn off Bluetooth when not in use and be cautious about connecting to unknown devices, especially in a public setting.

Range

The product must be within a range of approximately 33 feet of the device to which it is connected.

Expected Service Life

The expected service life of the product is three years.

If you dispose of the product, please dispose of it in accordance with your local, state, and/or country electronic recycling laws.

Intended Operators

Ordinary person, children likely present, instructed person, skilled person; teen or adult. Children under 12 may need adult supervision when connecting the speaker or controlling volume.

Expected Position of Operators

The product can be attached to a VersaWrap™. With the product attached, the wrap is expected to be within a variable distance of the individual's body.

Special Skills, Training, or Knowledge Required

No special skills or training are required to operate or maintain the product. The pictures and text in the instructions the manual will help the operator acquire the knowledge to achieve the intended use of the product.

Instructions for Use

- All functions of the product are safe for use by the user. For best results, follow the instructions provided with the product.
- For instruction on proper mounting of in a VersaWrap, see the instructions for use provided with the product.
- For information on cleaning the product, see "Cleaning, Maintaining, and Servicing" under "Safety Information".

Contact Information for Assistance

If you need assistance in setting up, using, or maintaining the product or if you need to report unexpected operation or an unexpected event related to the product, please contact versasupport@prc-salttillo.com or visit prc-salttillo.com/versa. For products purchased directly from Amazon, contact Amazon for support. To report any serious incident associated with the product, also contact your national authority.

Instructions for Storage

If you will be leaving the product in storage for a long period of time, be sure to turn all equipment off. Equipment with batteries in storage must be charged every six months.

Explanation of Symbols

This section explains all symbols that appear on all PRC-Salttillo products or packaging. To find an explanation for any symbol that appears on the product, scroll down through this table until you find it.

Symbol	Explanation
	The product is Federal Communications Commission (FCC)-compliant.
	The product conforms to European Union health, safety, and environmental standards.
	The product conforms to the requirements of the UK MDR 2002.
	The product complies with Innovation, Science and Economic Development Canada license-exempt RSS standard(s).
	The product complies with applicable ACMA technical standards for telecommunications, radiocommunications, EMC, and EME.
	Caution! Use caution when operating the product to avoid undesirable consequences.
	The product emits generally elevated, potentially hazardous, levels of non-ionizing electromagnetic radiation.
	Recycle electronic equipment. Do not throw the product in the trash.
IP22	The product is protected against solid objects 12.5 mm and greater and water dripping vertically at up to a 15-degree angle.
	Keep dry. The product needs to be protected from moisture.
	The range of temperatures to which the product can be exposed while in use or while in transport or storage. Caution: If the product has been stored in an environment in which the temperature range is different from the range of temperatures to which the product can be exposed while in use, let the product stand for at least 30 minutes before using it.
	The range of humidity to which the product can be exposed while in use.
	The range of atmospheric pressure to which the product can be safely exposed while in use, transport, or storage..
	Follow the instructions for use. Refer to the instruction manual.
	Consult the instructions for use.
	The product employs Bluetooth wireless technology.
	The product model number.
	The product serial number.
	The date of manufacture of the product and the name and address of the manufacturer.
	Indicates a medical product.
	The authorized representative in Switzerland.
	The authorized representative in the European Community.
	The responsible person in the United Kingdom.

International Contact Information



Manufacturer

Prentke Romich Company
1022 Heyl Rd
Wooster, OH 44691
USA
Tel.: +1 330 262-1984 · +1 800 262-1984
Fax: +1 330 262-4829
Email: info@prentrom.com
Web: www.prc-salttillo.com/



EU Authorised Representative

Prentke Romich GmbH
Karthäuserstr. 3
D-34117 Kassel
Germany
Tel.: +49 (0) 561 7 85 59-0
Fax: +49 (0) 561 7 85 59-29
Email: info@prentke-romich.de
Web : www.prentke-romich.de/



EU Importer

Prentke Romich GmbH
Karthäuserstr. 3
D-34117 Kassel
Germany
Tel.: +49 (0) 561 7 85 59-0
Fax: +49 (0) 561 7 85 59-29
Email: info@prentke-romich.de
Web : www.prentke-romich.de/



UK Responsible Person and Importer

Liberator Ltd
"Whitegates", 25-27 High Street
Swinstead, Lincolnshire NG33 4PA
UK
Tel.: +44 (0) 1733 370 470
Fax: +44 (0) 1476 552 473
Email: info@liberator.co.uk
Web: www.liberator.co.uk/

Authorised Representative in Australia

Liberator Pty Ltd
265, Gilbert Street,
Adelaide,
SA 5000
Australia
Tel: +61 (0) 8 8211 7766
Email: info@liberator.net.au
Web: www.liberator.net.au/