



VersaWrapTM

VersaWrap 8.1, VersaWrap 10.1, VersaWrap 12.1
Manual: English



DOC100322v101

Copyright and Trademark Information

© 2025 PRC-Salttillo. Non-commercial reprint rights for clinical or personal use granted with inclusion of copyright notice. Commercial use prohibited; may not be used for resale. Contact PRC-Salttillo for questions regarding permissible uses.

AAC Coaching Corner™, AAC Language Lab®, AAC Learning Journey™, Accent®, AeroSwitch®, AeroText™, ChatEditor™, ChatFusion™, ChatPoint™, ChatSMS®, ChatWrap®, CoreScanner™, Dialogue®, Empower®, Engage®, Essence®, Freedom Switch™, iSharePRC®, LAMP™, LAMP Words for Life®, Look®, Minspeak®, NovaChat®, NuEye®, NuPoint®, NuVoice®, PASS™, PRiO®, Realize Language™, Semantic Compaction®, Timocco-on-Accent™, TouchChat®, UNIDAD®, Unity®, VersaChat™, VersaEye™, VersaSpeaker™, VersaWrap™, Via®, VocabPC™, and Xtreme™ are trademarks of PRC-Salttillo.

AirDrop®, App Store®, Apple®, Apple Books®, FaceTime®, Guided Access®, iCloud®, iPad®, iPad Air®, iPad mini®, iPadOS®, iPad Pro®, iPhone®, iPod®, iTunes®, Lightning®, Mac®, Macintosh®, Safari®, Siri®, TouchID®, and TrueDepth® are trademarks of Apple Inc., registered in the U.S. and other countries and regions.

The *Bluetooth®* word mark and logos are registered trademarks owned by the Bluetooth SIG, Inc. and any use of such marks by PRC-Salttillo is under license.

iOS® is a trademark or registered trademark of Cisco in the U.S. and other countries and is used under license.

Wi-Fi® is a registered trademark of Wi-Fi Alliance.

All other trademarks and copyrights are the property of their respective owners.

Versa Support

Email: versasupport@prc-salttillo.com

Web: <https://www.prc-salttillo.com/versa>

Document Version

DOC100322v100 · 05/04/2026

Table of Contents

Safety Guidelines.....	3
VersaWrap Overview	4
Attach the VersaWrap to an iPad.....	5
Attach the Optional Carrying Strap.....	6
Attach the Optional VersaSpeaker.....	6
Specifications.....	7
Compliance Information	8
International Contact Information	10

Safety Guidelines

Cleaning, Maintaining, and Servicing

- ⚠ Do not attempt to disassemble, service, or repair the wrap. If the wrap is damaged, call Service for assistance.
- ⚠ Use a microfiber (lint-free) cleaning cloth to wipe the wrap. Never use abrasives or cleaning cloths that could damage it.
- ⚠ If the wrap needs more than a dry wipe-down, dampen a lint-free cloth. “Damp” means wrung out, never dripping. Dry the unit thoroughly with a dry lint-free cloth.

Water

- ⚠ This wrap is intended for use in normal situations. It is water-resistant, not waterproof. Do not immerse it in water. Do not spray or splash liquid directly on it and use extreme caution when using it around water. Never use the wrap when you are actually in the water (for example, a pool or bath).

Heat

- ⚠ Do not place the wrap in an appliance such as a microwave or oven, near an open flame, or on or near a source of heat such as a stove or heater.
- ⚠ Do not leave the wrap in a hot vehicle for a prolonged period.

Damage and Wear

- ⚠ Do not use the wrap if it is cracked or broken.
- ⚠ On a regular basis, visually inspect the wrap for wear or damage.

Environment

- ⚠ Dispose of all electronic components used with the wrap according to local regulations.

Attaching the VersaWrap™

- ⚠ Be sure to follow the instructions in this manual. Failure to assemble the wrap according to the instructions may result in an injury to the user.

Carrying Strap (sold separately)

- ⚠ The optional carrying strap is a potential entanglement or strangulation hazard. Please consider this prior to placing it with users.

VersaSpeaker (sold separately)

- ⚠ If you attach a VersaSpeaker™ to the VersaWrap, read and follow the Safety Guidelines in the VersaSpeaker Manual.

Apple AirTag (sold separately)


- ⚠ If you attach an Apple AirTag® to keep track of personal items from the VersaWrap, the AirTag poses security risks. Go to <https://support.apple.com/en-us/119874> and read the article “What to do if you get an alert that an AirTag, set of AirPods, Find My network accessory, or compatible Bluetooth location-tracking device is with you”.
- ⚠ If you attach an Apple AirTag to your VersaWrap, the AirTag, battery cover, and battery are potential choking hazards, especially to small children.

VersaWrap Overview

VersaWrap™ is a lightweight, free-standing wrap for iPad® made from an impact resistant material and featuring an easy-grip handle and multi-angle stand.


Options include a carrying strap, VersaSpeaker, and Apple AirTag®.

VersaWrap comes in three models: VersaWrap 8.1, VersaWrap 10.1, and VersaWrap 12.1. See “Specifications” on page 7 to determine the model you need for a specific iPad model.

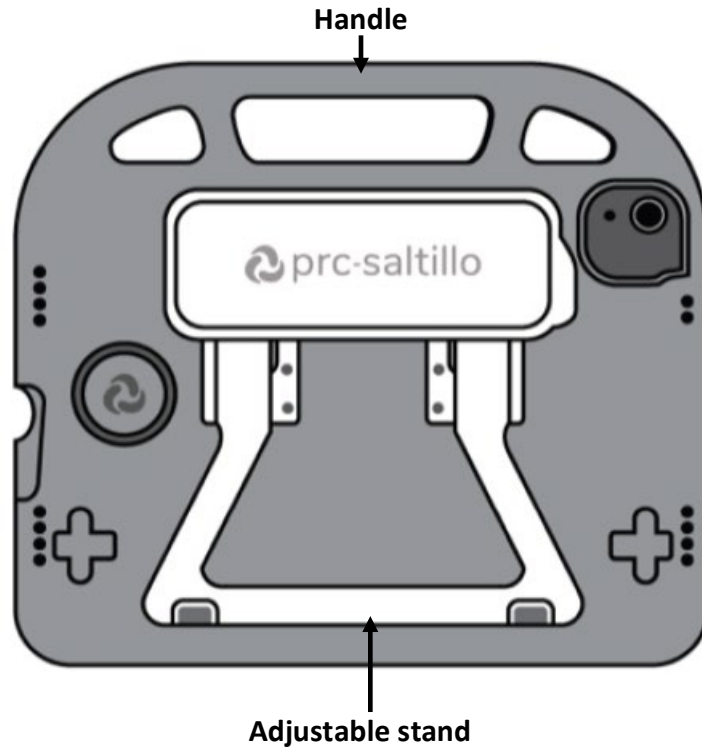


Attach the VersaWrap to an iPad

1. Remove the black plastic insert inside the VersaWrap.
2. Optional: Press an Apple AirTag (sold separately) into the built-in pocket.

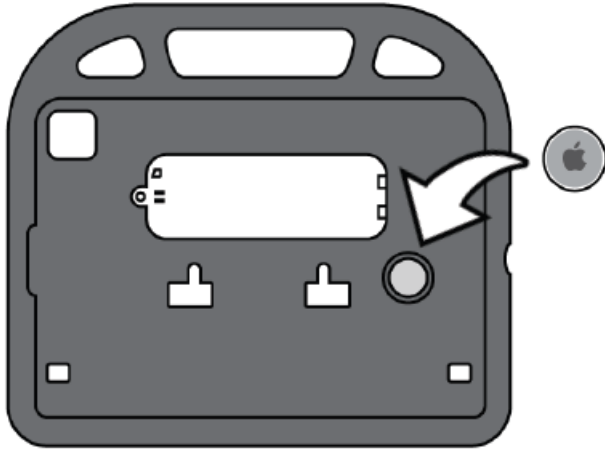


3. Insert the iPad toward the bottom right corner of the wrap.

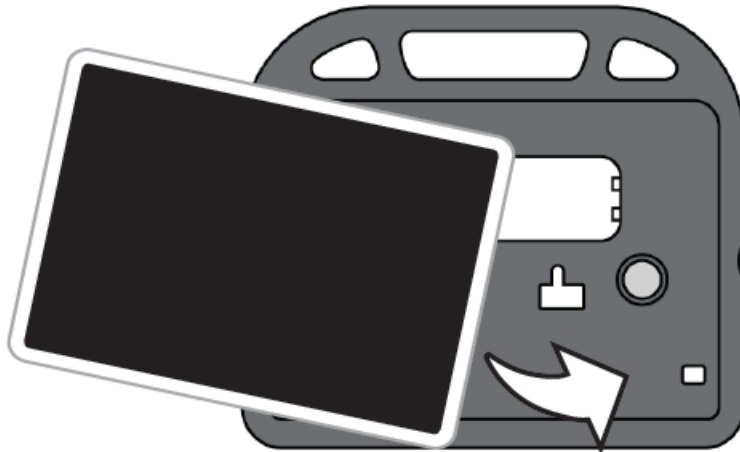

4. Carefully begin to pull the edges of the wrap onto the edges of the iPad until the iPad is fully enclosed by the wrap.


Attach the Optional Carrying Strap

If you want to carry the VersaWrap using a strap, attach the carrying strap through the eyelets on the back lower corners of the VersaWrap case.

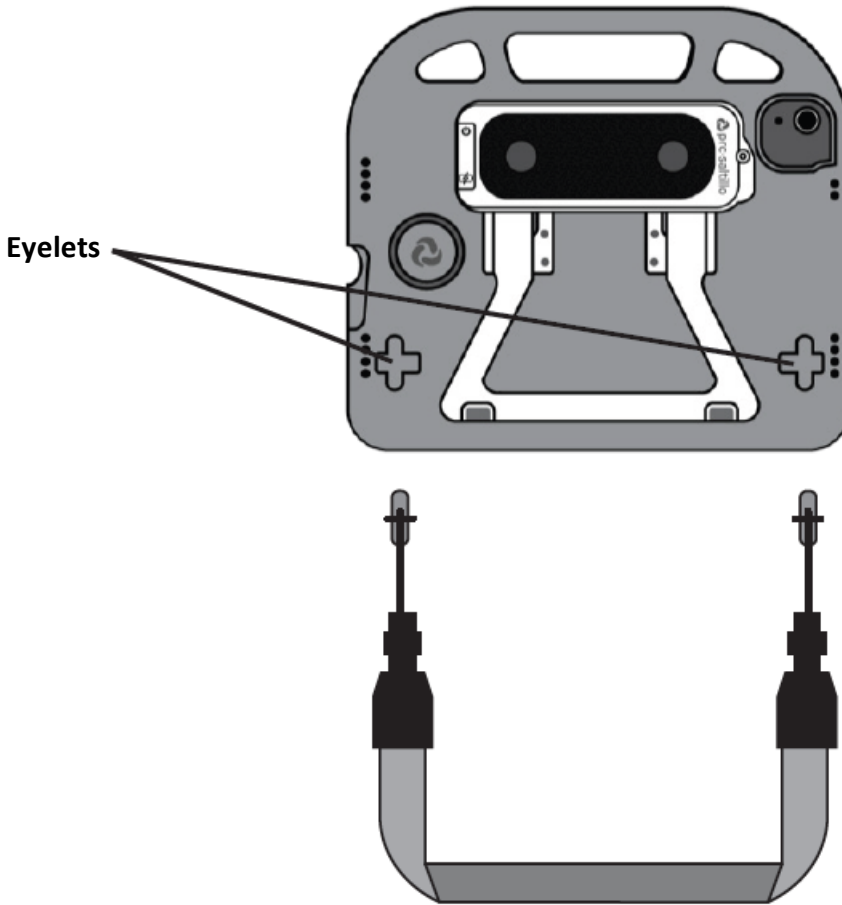

Attach the Optional VersaSpeaker

Amplify your sound with VersaSpeaker (sold separately). To attach a VersaSpeaker to the VersaWrap, follow the instructions in the Via Speaker User Guide.

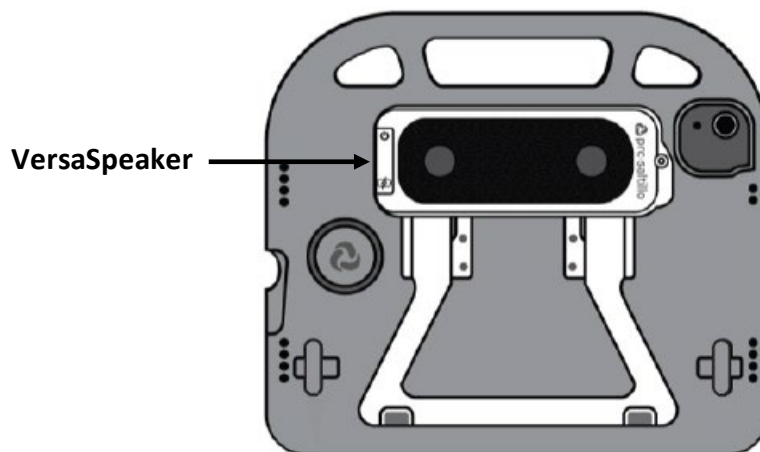

Specifications

VersaWrap 8.1

Weight

0.5 lbs.

Dimensions

7.5" L x 8.5" W x 1.7" D

Compatible with

iPad mini® 6th Gen and iPad mini (A17 Pro)

VersaWrap 10.1

Weight

0.6 lbs.

Dimensions

9.1" L x 10.4" W x 1.6" D

Compatible with

iPad Pro 11 1st, 2nd, 3rd, and 4th Gen and 11 (M4)

iPad Air 4th and 5th Gen and 11 (M2)

iPad 7th, 8th, 9th, and 10th Gen and 11 (A16)

VersaWrap 12.1

Weight

0.7 lbs.

Dimensions

10.6" L x 11.7" W x 1.7" D

Compatible with

iPad Pro 12.9 3rd, 4th, 5th, and 6th Gen and 13 (M4)

iPad Air 13 (M2)

Compliance Information

Product Intended Purpose

The product is intended to allow the user to protect, view, and carry an iPad. The optional carrying strap, VersaSpeaker, and Apple AirTag, provide additional functionality.

Contraindications:

- When an iPad is attached to the wrap, the fully assembled unit may be too heavy for some children under 12 to carry safely.
- Children under 12 may need adult supervision when connecting an optional carrying strap or VersaSpeaker.
- Attaching an optional Apple AirTag creates potential security issues.

Expected Service Life

The expected service life of the product is three years.

If you dispose of the product, please dispose of it in accordance with your local, state, and/or country electronic recycling laws.

Intended Operators

Ordinary person, children likely present, instructed person, skilled person; teen or adult. Some children under 12 may need adult supervision to use a fully assembled wrap safely.

Expected Position of Operators

The device can be viewed on a flat surface or carried with the handle or an optional carrying strap. A fully assembled unit is expected to be within variable distances of the individual's body.

Special Skills, Training, or Knowledge Required

No special skills or training are required to operate or maintain the product. The pictures and text in the instructions the manual will help the operator acquire the knowledge to achieve the intended use of the product.

Instructions for Use

- All functions of the product are safe for use by the user. For best results, follow the instructions provided with the product.
- For instructions on attaching an optional carrying strap, see the instructions provided with the strap.
- For instructions on proper mounting of a VersaSpeaker, see the instructions provided with that product.
- For information on cleaning the product, see "Cleaning, Maintaining, and Servicing" under "Safety Information".

Contact Information for Assistance

If you need assistance in setting up, using, or maintaining the product or if you need to report unexpected operation or an unexpected event related to the product, please contact versasupport@prc-salttillo.com or visit prc-salttillo.com/versa. For products purchased directly from Amazon, contact Amazon for support. To report any serious incident associated with the product, also contact your national authority.

Instructions for Storage

If you will be leaving the product in storage for a long period of time, be sure to turn all associated equipment off. Equipment with batteries in storage must be charged every six months.

Explanation of Symbols

This section explains all symbols that appear on all PRC-Salttillo products or packaging. To find an explanation for any symbol that appears on the product, scroll down through this table until you find it.

Symbol	Explanation
	The product is Federal Communications Commission (FCC)-compliant.
	The product conforms to European Union health, safety, and environmental standards.
	The product conforms to the requirements of the UK MDR 2002.
	The product complies with Innovation, Science and Economic Development Canada license-exempt RSS standard(s).
	The product complies with applicable ACMA technical standards for telecommunications, radiocommunications, EMC, and EME.
	Caution! Use caution when operating the product to avoid undesirable consequences.
	The product emits generally elevated, potentially hazardous, levels of non-ionizing electromagnetic radiation.
	Recycle electronic equipment. Do not throw the product in the trash.
	The product is protected against solid objects 12.5 mm and greater and water dripping vertically at up to a 15-degree angle.
	Keep dry. The product needs to be protected from moisture.
	The range of temperatures to which the product can be exposed while in use or while in transport or storage. Caution: If the product has been stored in an environment in which the temperature range is different from the range of temperatures to which the product can be exposed while in use, let the product stand for at least 30 minutes before using it.
	The range of humidity to which the product can be exposed while in use.
	The range of atmospheric pressure to which the product can be safely exposed while in use, transport, or storage..
	Follow the instructions for use. Refer to the instruction manual.
	Consult the instructions for use.
	The product employs Bluetooth wireless technology.
	The product model number.
	The product serial number.
	The date of manufacture of the product and the name and address of the manufacturer.
	Indicates a medical product.
	The authorized representative in Switzerland.
	The authorized representative in the European Community.
	The responsible person in the United Kingdom.

International Contact Information



Manufacturer

Prentke Romich Company
1022 Heyl Rd
Wooster, OH 44691
USA
Tel.: +1 330 262-1984 · +1 800 262-1984
Fax: +1 330 262-4829
Email: info@prentrom.com
Web: www.prc-salttillo.com/



EU Authorised Representative

Prentke Romich GmbH
Karthäuserstr. 3
D-34117 Kassel
Germany
Tel.: +49 (0) 561 7 85 59-0
Fax: +49 (0) 561 7 85 59-29
Email: info@prentke-romich.de
Web : www.prentke-romich.de/



EU Importer

Prentke Romich GmbH
Karthäuserstr. 3
D-34117 Kassel
Germany
Tel.: +49 (0) 561 7 85 59-0
Fax: +49 (0) 561 7 85 59-29
Email: info@prentke-romich.de
Web : www.prentke-romich.de/



UK Responsible Person and Importer

Liberator Ltd
"Whitegates", 25-27 High Street
Swinstead, Lincolnshire NG33 4PA
UK
Tel.: +44 (0) 1733 370 470
Fax: +44 (0) 1476 552 473
Email: info@liberator.co.uk
Web: www.liberator.co.uk/

Authorised Representative in Australia

Liberator Pty Ltd
265, Gilbert Street,
Adelaide,
SA 5000
Australia
Tel: +61 (0) 8 8211 7766
Email: info@liberator.net.au
Web: www.liberator.net.au/